# TenZ

---

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## 1. DEVICE DESCRIPTION

The Fergus-brand, **Tenz**-model, titanium button with adjustable loop comes as a metallic button with an adjustable loop, a blue #5 traction suture and a white-blue #2 flipping suture.

The adjustable loop is highly-resistant and is suitable for bone tunnels with different lengths, allowing for an easy graft placement. This product is indicated for soft-tissue-graft-to-cortical-bone suspension fixation.

## 2. PRODUCT PRESENTATION

Contents:

- One titanium button with adjustable loop, one blue #5 traction suture and one white-blue #2 flipping suture.

## 3. INDICATIONS FOR USE

The titanium button with adjustable loop is intended to offer suspension fixation of a graft to the cortical bone for anterior and posterior cruciate ligament reconstruction in the knee.

## 4. CONTRAINDICATIONS

- Insufficient bone quality or quantity for attachment. The physician must evaluate carefully the bone quality of patients with bone immaturity before starting with the orthopedic surgical procedure, and must avoid deviating, interrupting or altering the growth cartilage.
- Limited blood irrigation and/or previous infection which may delay consolidation.
- Conditions that may limit the ability and/or the will of the patient to follow the physicians' indications regarding post-operative care, which may include significant restrictions on their physical activity.
- Inappropriate neuromuscular conditions.

## 5. WARNINGS

- Read the instructions for use carefully and completely before use.
- The device must be used only by a medical professional trained in the surgical techniques indicated for placing these implants.
- The law restricts the use of this device to physicians or people authorized by a physician.
- DO NOT use the device if the packaging is open or damaged since implants are provided sterile.
- The device must not be used for surgical procedures that are not included in the “Indications for Use” section.
- Internal fixation devices have been designed for SINGLE use only. Therefore, THEY MUST NOT BE REUSED, REPROCESSED, REASSEMBLED OR RESTERILIZED, since this could impair their efficiency and increase the risk of inappropriate sterilization and cross-contamination.
- This is a single-use device. Reusing the device could lead to improper functioning and harm the patient and/or user.
- Do not use additional sutures. Any additional sutures could prevent the device from passing through the bone tunnel.
- The strictest precautions must be implemented to prevent contamination. Operating room conditions must be according to local hospital, administrative or governmental requirements.
- Do not use the device after the expiration date that appears in the label. Device performance, safety and/or sterility cannot be guaranteed after the expiration date has passed.
- The product must be stored in the original packaging.
- Patient sensitivity to device materials must be considered before implantation.
- All surgical instruments used in the surgical procedure to install this implant must be sterile and suitable for use in aseptic environments.
- After the procedure and until the complete bone consolidation is achieved, the fixation provided by the device must be considered as temporary. Therefore, it may not be able to resist weight loads or other efforts without assistance. The fixation provided by the device must be protected. Post-operative guidelines prescribed by the physician must be followed carefully in order to prevent overstressing the device and the bone.

- If the device is to be explanted, the risks of the possibility of a second procedure on the patient must be considered. After removing the implant, appropriate post-operative controls must be carried out.
- Make sure that the product is disposed of as per local regulations, and be wary of the risks of environmental, staff and patient contamination.
- Bio-hazardous waste, such as explanted devices, needles, sharp objects and contaminated surgical materials, must be disposed of safely and as per the policies of the site where the surgical procedure is performed.
- The manufacturer shall not be held responsible for faults or damages caused when these devices are combined with products not authorized/recommended by the company.
- The manufacturer shall not be held responsible for any effects or consequences that may arise from deviating from the specific surgical techniques or instructions, which may include, but not limited to, an incorrect indication or a misuse of the surgical technique, or problems with asepsis.
- The manufacturer does not recommend the use of Kocher forceps or strong (or traumatic) grasping forceps to hold the sutures when raising the graft or the implant. Using such forceps may damage the suture fibers and could even cut them.

## 6. PRECAUTIONS

- Surgeons are advised to review the product-specific surgical technique before the procedure.
- When material sensitivity is suspected, the corresponding tests must be carried out and sensitivity must be ruled out before implantation.
- Before using, visually inspect the device for damage. Do not use if damaged.
- Sutures may become damaged if handled or controlled using sharp or cutting instruments.
- Pre-operative and post-operative procedures, including the knowledge about surgical techniques and proper device selection and installation, are essential for satisfactory use of the device.
- Knee position must be kept immobilized. Knee flexion angle must not change during device insertion.
- Procedures before the operation and during the operation, including the knowledge of the surgical technique and

proper selection and positioning of the device, are fundamental for satisfactory implant use.

## **7. SUPPLY AND STORAGE**

- The titanium button with adjustable loop is provided sterile and pyrogen-free. It is sterilized using ethylene oxide. In order to maintain product quality and effectiveness, it is marketed in a double Tyvek pouch. The product must be used only if the original packaging and labeling are not damaged.
- The device must be kept in its original packaging, closed and away from any source of humidity.
- Room temperature storage is advised for this product.
- **IF THE PACKAGING IS OPEN OR ALTERED, DO NOT USE THE PRODUCT.**
- **DO NOT USE THE PRODUCT AFTER THE EXPIRATION DATE PRINTED IN THE PACKAGING.**
- The implant has a 3-year shelf life.

## **8. MAGNETIC RESONANCE (MR) ENVIRONMENTS**

### **1. Conditional MR Compatibility**

The titanium button with adjustable loop is made with the Ti-6Al-4V (ASTM F136) titanium alloy, which has a low magnetic susceptibility and electrical conductivity. This reduces the risk of displacement caused by magnetically-induced force or torque, radiofrequency-induced heating, and imaging artifacts.

Patients with an implanted titanium button with adjustable loop may be subjected to magnetic resonance screenings using the following conditions:

- Static magnetic field of 1.5 tesla and 3 teslas, only.
- Maximum spatial gradient magnetic field of 3000 gauss/cm or less.
- Maximum whole-body averaged specific absorption rate (SAR) reported by the MR system of 2 W/kg for 15 minutes of scanning in the normal mode of operation for the MR system.

Image quality could be affected in the implant area or its surroundings. Therefore, magnetic resonance parameters may need to be optimized to compensate for the presence of the device.

## **2. MR Safety**

Fergus-brand sutures, manufactured with UHMWPE, are MR-safe.

## **9. POSSIBLE COMPLICATIONS**

- Infection, both deep and superficial.
- Pain on the site of the incision or the procedure.
- Inflammatory or allergic reactions.
- Reactions to foreign bodies.
- Implant loosening and/or detachment.
- Bone lesions.
- Incomplete or inappropriate correction, or lesion recurrence.
- After surgery, ligament laxness can be affected, and therefore, the knee range of motion may decrease or increase.

## **10. PATIENT INFORMATION**

The surgeon is responsible for informing the patient and his/her representatives before surgery about the possible complications related to the implantation of the device. The manufacturer and its distributors delegate to the surgeon the responsibility of informing the patient about the advantages and possible risks associated with the implantation and use of the device.

Post-operative care is important. The patient will be instructed about the limitations required due to the implant and will be warned about weight unloading and device overloading before a safe bone consolidation.

The patient must be warned that, if explantation of the device is deemed necessary, the risks of the possibility of a second surgical procedure on the patient must be considered. Device extraction must be accompanied by an appropriate post-operative approach.

Performing additional surgeries does not guarantee a complete repair of the complications.

## **11. SURGICAL TECHNIQUE**

The Tenz titanium button with adjustable loop is used in the anterior or posterior cruciate ligament reconstruction for the fixation of soft tissue graft to a femoral or tibial bone

tunnel, as per surgeon decision.

The following surgical technique instructions refer only to device implantation, not considering other procedures that are critical to perform the surgical procedure, such as patient preparation, working portals, procedure finalization, etc.

Bone tunnels may be performed by using the following surgical techniques:

- Transportal technique (through the anteromedial arthroscopic portal)
- Transtibial technique (through a tibial bone tunnel)
- All-inside technique

General steps:

**1.** Remove and condition the graft by using the appropriate techniques and instruments.

**2.** Mount the graft over the adjustable loop of the implant on the graft preparation table.

Prepare the remaining end of the graft as required for the preferred fixation method.

**3.** Create the bone tunnels and channels considering the chosen surgical technique, the anatomical references and the dimensions of the graft, and using the appropriate instruments.

**A passage channel with a minimum diameter of 4 mm to a maximum diameter of 5 mm is needed to pass the button through the bone and achieve optimal and safe fixation.**

When creating the bone tunnel where the graft is to be housed, make sure that a bone channel with the appropriate length is maintained. The length should not be less than 10 mm.

**Note:**

- Mark the implant at a distance equal to the total length of the condyle, measuring from the metallic button towards the loop. This marking is used when raising the button, in order to determine if it is already in the cortical bone.
- Mark the graft at a distance equal to the length of the femoral tunnel, measuring from the beginning of the intersection between the graft and the Tenz loop. This marking is used when raising the graft in order to determine if it is taking up the entire length of the bone tunnel.

**4.** Pass the shortening suture, the blue #5 traction suture and the white-blue #2 flipping suture through the bone channel. Place the graft inside the joint using the bone tunnel created for the task. Pull the traction suture of the

device until the button goes through the bone channel completely in order to pass the graft through the tunnels, and use the flipping suture to position the button so that it lays completely on the cortical surface of the bone.

**5.** Once the button is placed correctly on the cortical bone surface, pull both ends of the adjustable loop, while maintaining tension on the distal end of the graft, until the desired tension is achieved.

Note: The ends of the adjustable loop must be symmetrically tensioned.

**6.** Remove the traction and flipping sutures, and cut the excess length of the adjustable loop.

Note: Avoid cutting the shortening sutures next to the splice in the button, since this may compromise fixation.

Precautions:

- Graft diameter must be measured accurately to make sure it matches the diameter of the bone tunnel and avoid complications with the bone tunnel, the suture and/or the graft during the surgical procedure.
- Applying tension to the white-blue #2 flipping suture while the button is being passed through the bone tunnel could flip the device to a position perpendicular to the graft, creating unwanted friction on the bone tunnel.
- No additional force needs to be applied to the shortening suture strands when the graft construction reaches the desired position in the femoral cavity and graft stability is verified by pulling from it distally.

## 12. MATERIAL SPECIFICATIONS

The titanium button with adjustable loop is manufactured using the following materials:

- The button is made with a titanium alloy, Ti-6Al-4V ELI, as per ASTM F136.
- The sutures are manufactured with ultra-high-molecular-weight polyethylene (UHMWPE).

Check the packaging labels for information about the materials.


## 13. LEARN MORE

To learn more about this product, visit our website:  
[www.fergusmedical.com](http://www.fergusmedical.com)

## 14. SYMBOLS USED IN LABELS

**CATALOG NUMBER**

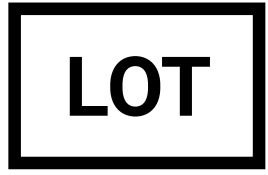
**BATCH NUMBER**


**SERIAL NUMBER**


**PRECAUTION**

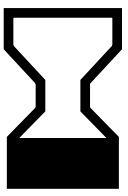**EXPIRATION DATE**

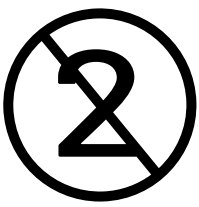**DO NOT REUSE**

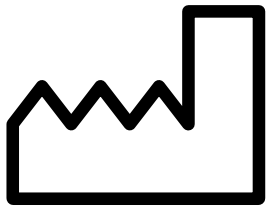
**MANUFACTURING DATE**


**STERILE STERILIZATION METHOD: ETHYLENE OXIDE**

**MANUFACTURER**

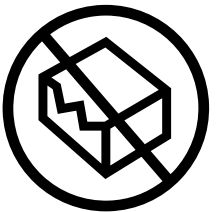
**DO NOT USE IF PACKAGING IS DAMAGED**

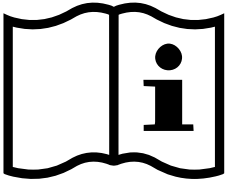**CHECK INSTRUCTIONS FOR USE**

**DO NOT RE-STERILIZE**

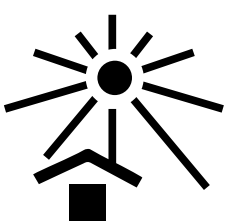**KEEP AWAY FROM SUNLIGHT**

**KEEP DRY**

**MR-CONDITIONAL**

# Fergus

**PROMEDON S.A**

Av. Gral. Manuel Savio s/n

Lote 3 • Manzana 3 • (X5123XAD)

Pque Ind. Ferreyra, Córdoba, Argentina

[www.fergusmedical.com](http://www.fergusmedical.com)